



ReSound GN

ReSound ONE

Hear like no other

GN Making Life Sound Better

FOR PROFESSIONALS

The world's first full-featured* hearing aid with M&RIE

JUST LIKE A FINGERPRINT, EACH PERSON'S HEARING IS ONE OF A KIND

Nature has designed the shape of the ear to perfectly capture the sounds around it. And because the shape of the ear is unique to the individual, the sound they hear is uniquely theirs. Unfortunately, when wearing hearing aids, this isn't the case, as most devices capture sound from behind the ear only. Needless to say, this isn't where nature intended, so sound ends up feeling standardised, unnatural and incomplete.

ReSound ONE™ with M&RIE (Microphone and Receiver-In-Ear) changes all of this, an entirely new class of hearing aids* with a microphone and receiver positioned inside the ear, so users can collect the whole sound, just as nature intended. This breakthrough is made possible by a brand-new sound processing chipset and feedback management system. ReSound ONE gives the brain everything it needs to process sound with greater depth and direction, making it effortless to tune in and enjoy conversation, in any situation.

This is the next step to what we call Organic Hearing.

ReSound ONE. Hear like no other.

*Microphone in ear and two standard directional microphones, directionality features, wireless audio streaming



The ReSound Organic Hearing philosophy


Since its beginning, the ReSound sound philosophy has always been to emulate the natural hearing process as close as possible. We call that **Organic Hearing**. It's the foundation of our technology that leverages the individual's hearing anatomy to collect and direct sound organically to provide the brain with the most complete sound experience – even in the most challenging listening situations.

Everything you need to give your client the most individualised hearing experience



OUR MOST NATURAL SOUND QUALITY

The groundbreaking M&RIE* (Microphone & Receiver-In-Ear), inspired by our Organic Hearing philosophy, combines the two traditional microphones with a third that sits discreetly inside the ear canal as an integrated part of the receiver module.

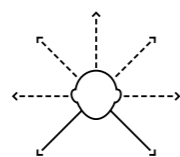

ALL ACCESS DIRECTIONALITY



Automatically adjusts directional microphone patterns to put your clients in the best position to hear what's important and still monitor the sounds around them.



ULTRA FOCUS

The user-controlled Ultra Focus uses a powerful directional pattern to focus on the person in front of your clients in very noisy environments.

*Microphone in ear and two standard directional microphones, directionality features, wireless audio streaming

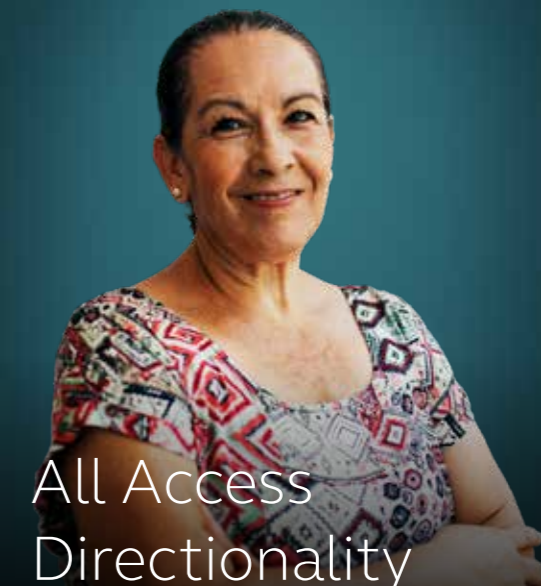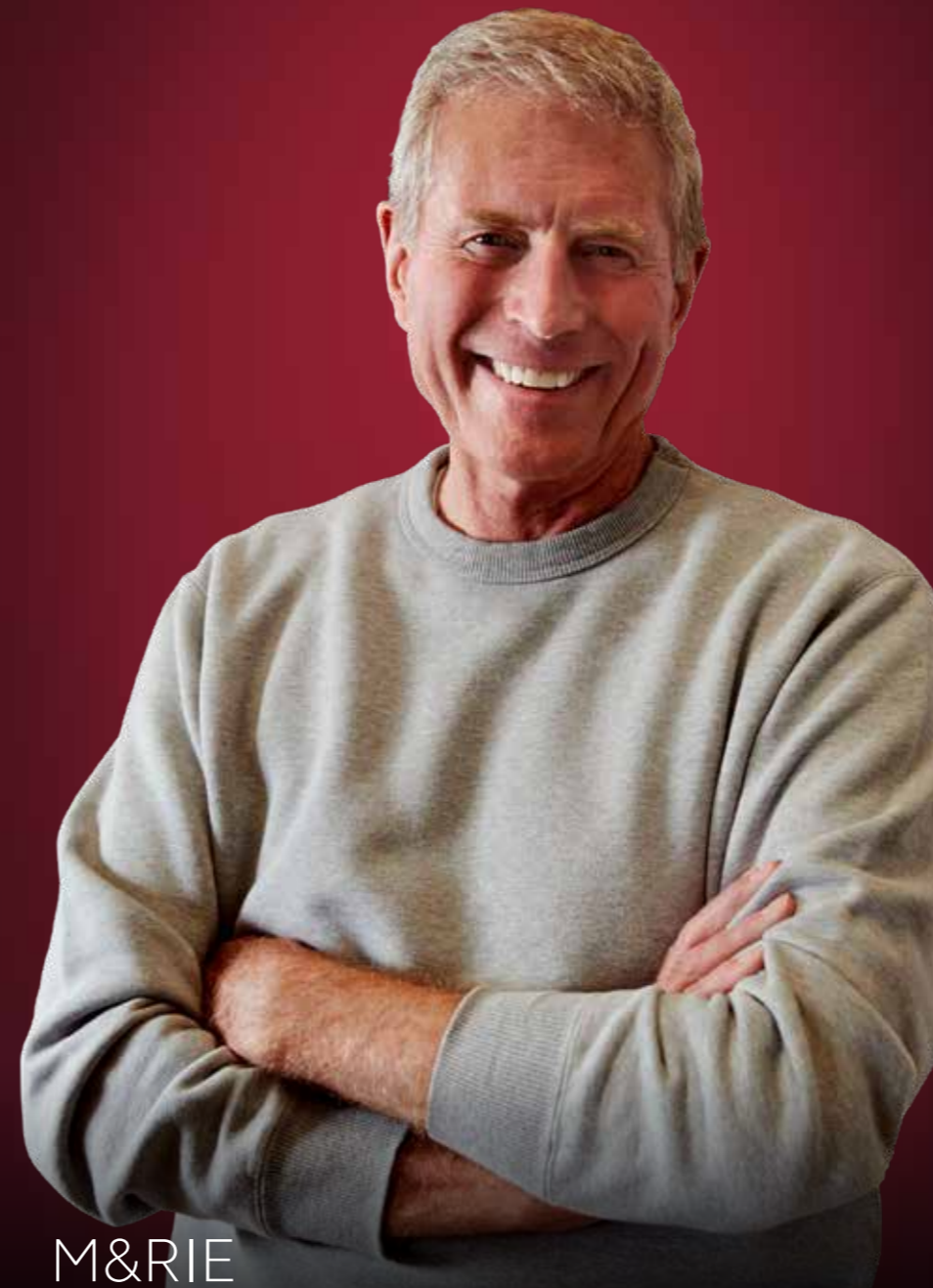



Our most natural sound quality in an RIE hearing aid

The M&RIE combines the two traditional microphones with a third in the ear canal to help clients hear with more direction and depth.


BEST OF BOTH WORLDS

Made possible by a new, faster sound processing chipset and a new feedback management system, DFS Ultra III, this breakthrough uses the unique shape of the ear to collect sound the way nature intended.

Clients get the natural sound of custom products, combined with the advanced features and ease of use of a Receiver-In-Ear model. In fact, 90% of listeners prefer M&RIE for natural sound quality, and M&RIE users recognise sounds from the front and back 10% more accurately.

REDUCED WIND NOISE

Having the microphone sit inside the ear canal also creates natural wind noise reduction by 15dB without reducing gain, so whether they are out for a run or just enjoying a summer day, clients can hear comfortably.

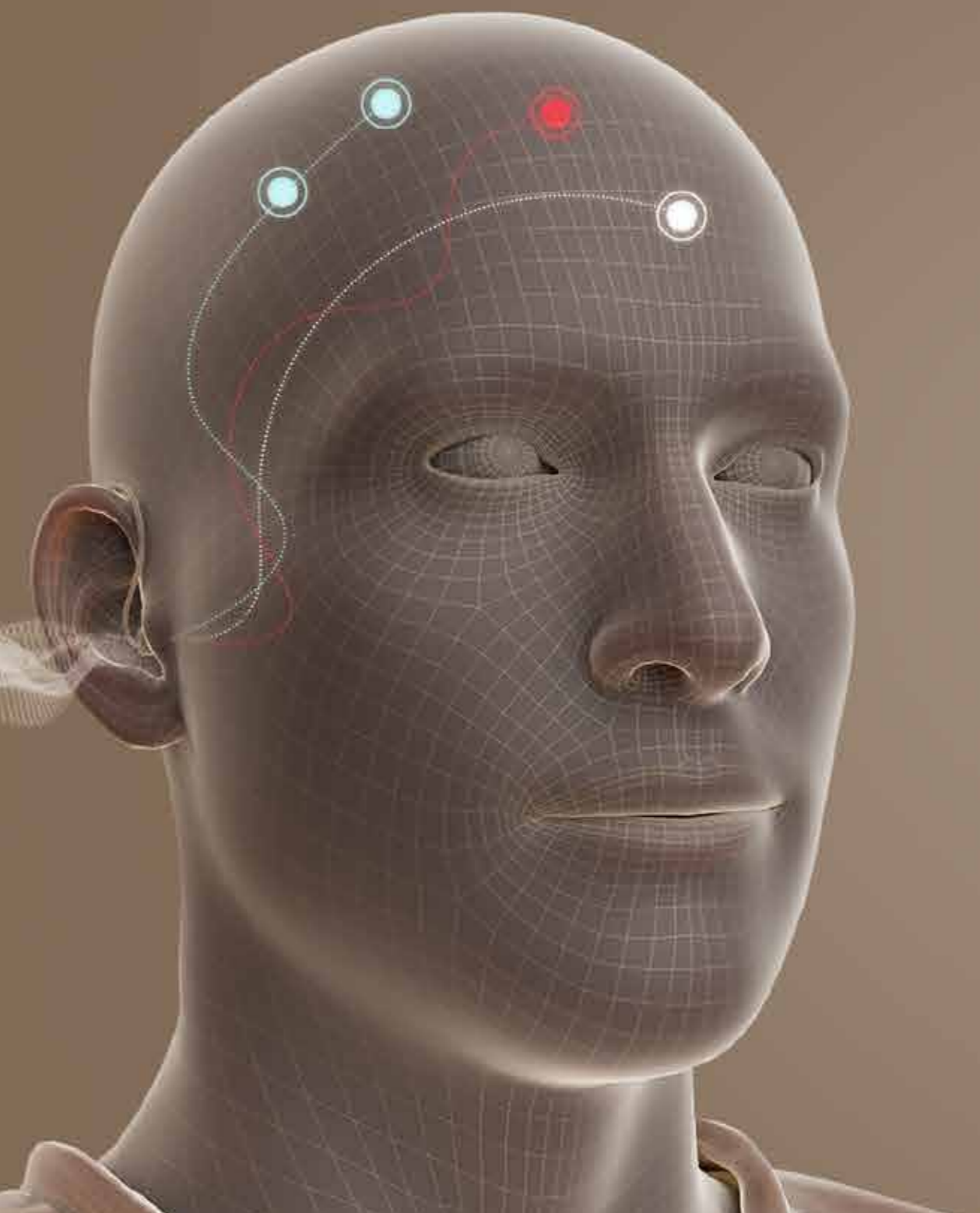


90% of listeners prefer* M&RIE for natural sound quality.

10% more accurate sound** recognition from the front and back.

15dB natural wind noise*** reduction without reducing gain.

*With mild hearing loss, and compared to pinna compensation and omnidirectionality
**With mild hearing loss, compared to traditional microphone configurations
***Compared to omnidirectionality
Source: Groth J. (2020). An innovative RIE with microphone in the ear lets users "hear with their own ears". ReSound white paper.

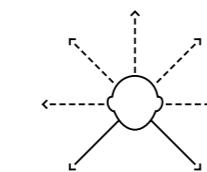

All Access Directionality

Greater hearing in any environment

Inspired by our Organic Hearing philosophy, All Access Directionality automatically adjusts directional microphone patterns to put your clients in the best and most natural position to hear what's important and still monitor the sounds around them. By allowing them to access more sounds than most hearing aids, it lets their brains, not their hearing aids, classify what is important, so they can feel at ease and engage effortlessly with others and their changing surroundings.

BETTER PERFORMANCE IN NOISE

In noisier situations, All Access Directionality utilises a new binaural connection to create a stronger, more focused beam for sounds in front of the listener. Compared to traditional directionality, it provides up to 2dB better signal-to-noise ratio.

ALL ACCESS DIRECTIONALITY AND M&RIE WORKING AS ONE

ReSound ONE monitors the listening situation and automatically switches among M&RIE and the different directional modes in All Access Directionality to ensure that your clients have the most optimised settings to hear their best in any situation.

2dB

better signal-to-noise ratio,
compared to traditional
directionality.*

*Vs. our legacy products, depending on hearing loss
Source: Groth J. (2020). The evolution of the binaural hearing strategy: All Access Directionality and Ultra Focus. ReSound white paper



Ultra Focus

Great one-to-one hearing experience in the toughest environments

The user-controlled Ultra Focus uses a powerful combined directional pattern from both hearing aids to focus on the person in front of your client. When they need extra help in very noisy environments, it provides a 30% improvement in speech understanding over All Access Directionality. Ultra Focus prioritises sound from the side with the least noise to make the speech as clear as possible. This unique feature gives up to 1.5dB signal-to-noise ratio improvement when noise is louder on one side. While speech is made clearer, Ultra Focus can still help clients maintain spatial cues around them, by using the speech frequencies for the targeted directionality in front and the high and low frequencies to monitor the environment.

EASY TO ACTIVATE WHEN THEY NEED IT MOST

Clients can choose Ultra Focus via the ReSound Smart 3D™ app or the program button on their ReSound ONE hearing aids for strong directionality. It's ideal for one-to-one conversations in very noisy situations, like a food hall, taxi ride with music and traffic noise, or service counter at a crowded airport. Ultra Focus can give clients renewed confidence in situations, which they might otherwise have avoided.


*Source: Groth J. (2020). The evolution of the binaural hearing strategy: All Access Directionality and Ultra Focus. ReSound white paper.

30%

improvement in speech understanding over All Access Directionality.*

1.5dB

signal-to-noise ratio improvement when noise is louder on one side.*

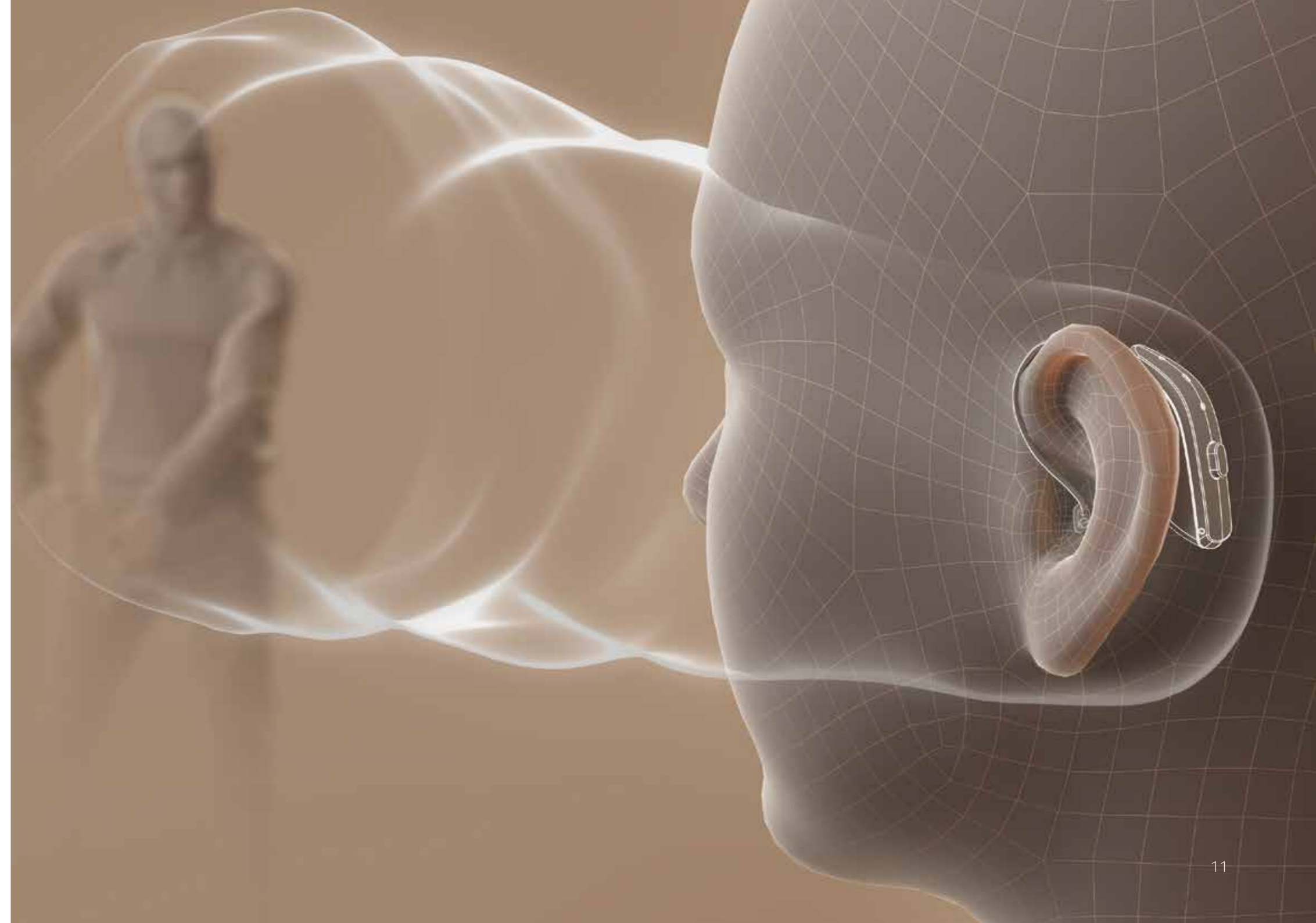



RIE62 13 Battery model

RIE61 Rechargeable model

RIE61 312 Battery model

Adapts to your clients' needs

With a new modern design and new range of attractive colours, new receivers with enhanced dome security, wireless calls, audio streaming from iPhone and Android™ smartphones, and great rechargeable battery life for worry-free all day usage, plus a more accurate, comfortable, and individualised first fit, ReSound ONE adapts to your clients' lives and makes it easy for them to connect to the world around them.



Warm Grey

Sparkling Silver

Graphite



Deep Black

Designed with a purpose - your client



Gold




Espresso



Champagne



Bronze

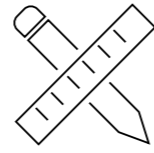

Design

Designed for ultimate comfort, discretion and performance

COMFORT AND DISCRETION

Every detail of ReSound ONE was designed with one purpose – the unique individual. Built for ultimate adaptability, the modern styling has just one comfortable touchpoint on the back of the ear, sits lower for greater discretion, and comes in a variety of colours to blend with hair and skin tones.


GREAT PERFORMANCE

For added flexibility, two different charger options serve as stylish, portable and protective cases, ideal for travelling and storage. The Premium Charger includes an onboard battery for up to three days of on-the-go portable charging, without the need for a power outlet, and the smaller Standard Charger gives clients the power they need in a simple, timeless design.

The new receiver design with enhanced dome security makes it easier to fit and easier for your clients to wear. Designed to the IP68 standard, with iSolate nanotech coating inside and out, ReSound ONE is resistant to water, dust and dirt, so they can rely on it throughout every day.



Premium Charger

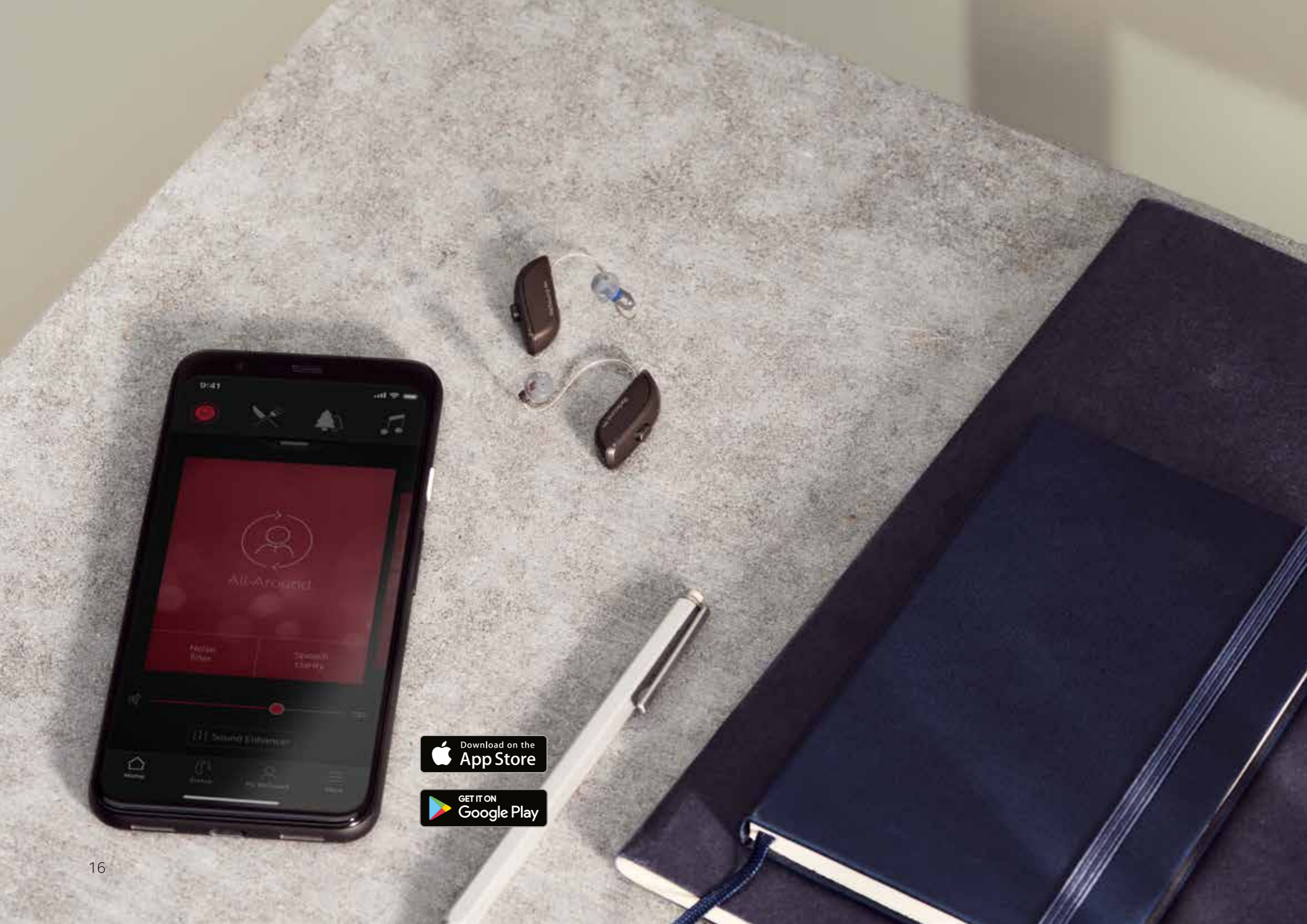

Download on the
App Store

GET IT ON
Google Play



Rechargeability & Connectivity

Unmatched rechargeability and extensive connectivity to adapt to clients' needs


THE POWER THEY NEED

ReSound ONE offers a rechargeable solution with the longest battery life in the industry, for well over a day of worry-free use. The 312 and 13 battery versions offer improved performance over previous models, so no matter which one you choose for your clients, they can be assured that ReSound ONE keeps working for them throughout their entire day.


A COMPLETE WIRELESS ECOSYSTEM

The long battery life enables ReSound ONE to adapt to clients' lifestyles as it simply connects them to their iPhone and Android smartphones and the world around them with wireless calls and direct audio streaming, via the stable Bluetooth® Low Energy.

Access to the ReSound Smart 3D app and a versatile range of wireless accessories give them even more options to connect to the sounds they enjoy every day.

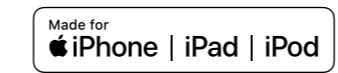

Phone Clip+ Multi Mic Micro Mic TV Streamer 2 Remote Control 2 Remote Control

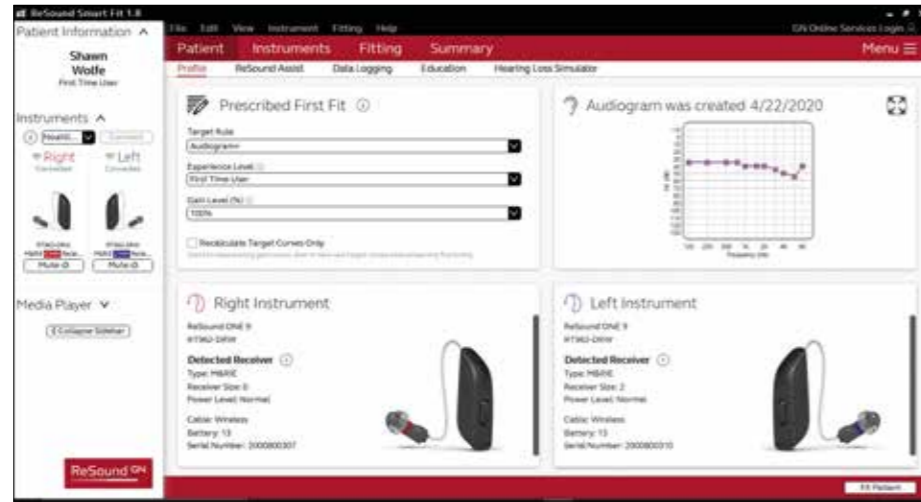


ReSound Smart Fit

An accurate, comfortable, and individualised first fit and beyond

The new and unique personalised feature prescriptions in ReSound Smart Fit™ are guided by peer-reviewed research on how to optimise advanced sound processing features. It can give you a better starting point ReSound ONE, to help decrease the need for fine-tuning. The ability to adjust gain percentage also helps you create the most comfortable first fit for your clients. A better start can lead to more client satisfaction and less returns.

EASE THROUGH FITTINGS

The ReSound Smart Fit flow is more intuitive than ever with Automatic Receiver Detection. When you connect ReSound ONE, the receiver is detected and the receiver specific calibration data is applied to the fitting. ReSound Smart Fit guides you through the process to make fittings more efficient.

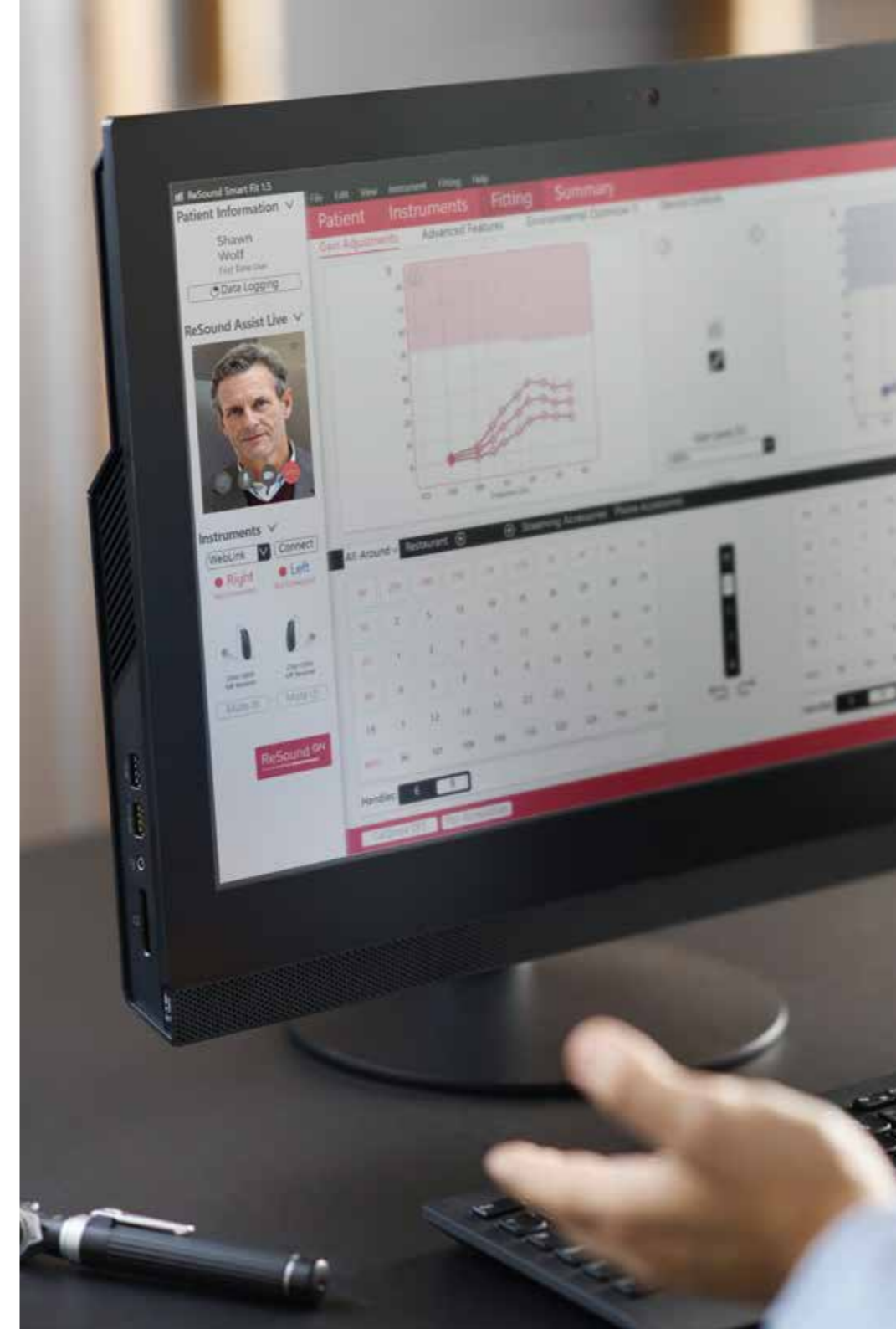

ReSound Assist and ReSound Assist Live

More connections for better care

With ReSound Assist and ReSound Assist Live, ReSound ONE gives you two powerful ways to serve and engage with your clients to provide better care, even if circumstances do not allow clients to be onsite with you. You'll be able to deliver more accessibility and the best possible fitting, counselling and support for your clients in the most convenient and safest way possible.

WHEN CLIENTS CAN'T GET TO YOU, BRING THE CARE TO THEM

When clients are unable to travel to your clinic, ReSound Assist Live video appointments give you a convenient solution for remote fine-tuning and fitting via ReSound Smart Fit. It's as easy as answering a call on their iPhone or Android smartphones. It includes at-home in-ear testing capabilities to give you the ability to offer the broadest range of services. You can send a trial set of hearing aids directly to their home, and with remote support, you can conduct a hearing assessment and immediately fit and fine-tune new hearing aids, ensuring the best possible individualised sound. When possible, clients can follow up with an in-clinic consultation.

Model M&RIE (MM Receiver)	RT961-DRW RT961-DRWC RT962-DRW	RT761-DRW RT761-DRWC RT762-DRW
Device Configurations		
Battery size 61-DRW	312 Zinc-Air	
Battery size 61-DRWC	Rechargeable Lithium-Ion	
Battery size 62-DRW	13 Zinc-Air	
Receiver Options	MM	
Control Options	PB (61-DRW and 61-DRWC), Multi-Function button (62), Telecoil (62)	
IP Classification	IP68	
Audiological Features		
Number of Channels	17	14
All Access Directionality	●	-
M&RIE	●	●
Ultra Focus	●	-
Synchronised Soft Switching	●	●
Environmental Optimiser II	●	-
Environmental Optimiser I	-	●
Noise Tracker II	5 settings	3 settings
Expansion	3 settings	2 settings
Impulse Noise Reduction	3 settings	2 settings
Wind Guard	3 settings	2 settings
DFS Ultra III (w/ Music Mode)	●	●
Synchronised Acceptance Manager	●	●
Tinnitus Sound Generator	●	●
Functional Features		
Ear to Ear Communication	●	●
Direct audio streaming	●	●
ReSound TV Streamer 2, Remote Control, Remote Control 2, Phone Clip+, Micro Mic and Multi Mic	●	●
ReSound Smart 3D app	●	●
Sound Enhancer (ReSound Smart 3D app)	●	-
ReSound Assist		
Remote Fine Tuning	●	●
Remote Firmware Updates	●	●
ReSound Assist Live	●	●
Fitting Features		
Fitting Software ReSound Smart Fit 1.8 or higher	●	●
Max Gain Handles	17	14
Fully Flexible Programs	4	4
Datalogging	●	●
Wireless Fitting with Noahlink Wireless	●	●
Automatic Receiver Detection	●	●



RT61-DRW



RT61-DRWC



RT62-DRW

Model (Standard Receivers)	RT961-DRW RT961-DRWC RT962-DRW	RT761-DRW RT761-DRWC RT762-DRW	RT561-DRW RT561-DRWC RT562-DRW	RT461-DRW RT461-DRWC RT462-DRW
Device Configurations				
Battery size 61-DRW	312 Zinc-Air			
Battery size 61-DRWC	Rechargeable Lithium-Ion			
Battery size 62-DRW	13 Zinc-Air			
Receiver Options	LP, MP, HP, UP			
Control Options	PB (61-DRW and 61-DRWC), Multi-Function button (62), Telecoil (62)			
IP Classification	IP68			
Audiological Features				
Number of Channels	17	14	12	12
All Access Directionality	●	-	-	-
Spatial Sense	●	●	-	-
Ultra Focus	●	-	-	-
Binaural Directionality III	-	●	-	-
Binaural Directionality	-	-	●	-
Synchronised Soft Switching	●	●	●	●
Environmental Optimiser II	●	-	-	-
Environmental Optimiser I	-	●	●	-
Noise Tracker II	5 settings	3 settings	2 settings	On/Off
Expansion	3 settings	2 settings	On/Off	On/Off
Impulse Noise Reduction	3 settings	2 settings	On/Off	On/Off
Wind Guard	3 settings	2 settings	On/Off	On/Off
Sound Shaper	●	●	●	●
DFS Ultra III (w/ Music Mode)	●	●	●	●
Synchronised Acceptance Manager	●	●	●	●
Tinnitus Sound Generator	●	●	●	●
Functional Features				
Ear to Ear Communication	●	●	●	●
Direct audio streaming	●	●	●	●
ReSound TV Streamer 2, Remote Control, Remote Control 2, Phone Clip+, Micro Mic and Multi Mic	●	●	●	●
ReSound Smart 3D app	●	●	●	●
Sound Enhancer (ReSound Smart 3D app)	●	-	-	-
ReSound Assist				
Remote Fine Tuning	●	●	●	●
Remote Firmware Updates	●	●	●	●
ReSound Assist Live	●	●	●	●
Fitting Features				
Fitting Software ReSound Smart Fit 1.8 or higher	●	●	●	●
Max Gain Handles	17	14	12	10
Fully Flexible Programs	4	4	4	4
Auto DFS	●	●	●	●
Datalogging	●	●	●	●
Wireless Fitting with Noahlink Wireless	●	●	●	●
Automatic Receiver Detection	●	●	●	●



RT61-DRW






RT61-DRWC



RT62-DRW

Fitting range

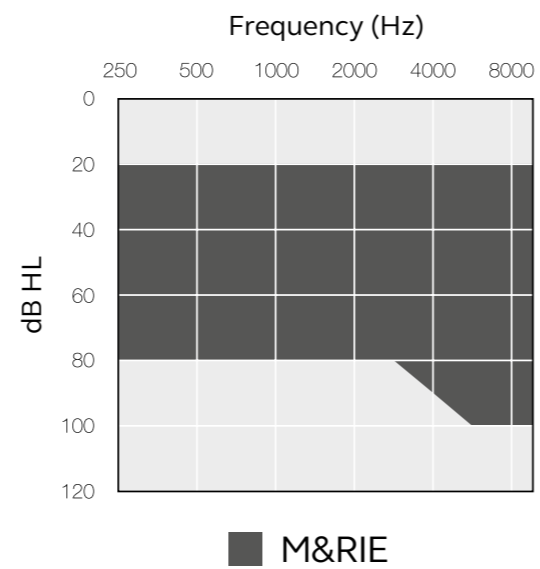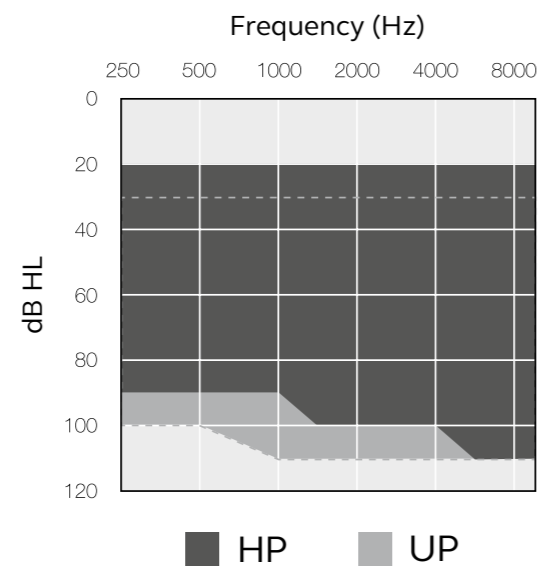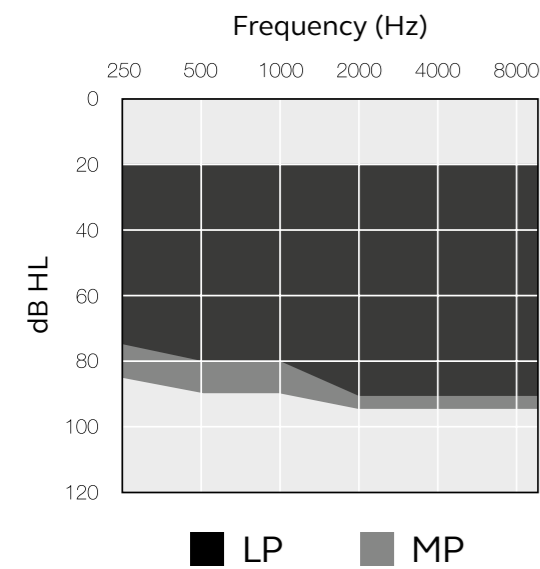

Chargers




Technical data	Premium Charger	Standard Charger
Dimensions	99.4 x 35 x 67.5 mm / 3,9 x 1.4 x 2.7"	100.2 mm x 42 mm x 54.8 mm / 3.9 x 1.7 x 2.2"
Weight	145 gram / 5.1 oz	95 gram/ 3.3 oz
Power Supply	USB power supply, 5 V	N/A
Power Connector	USB	USB
Internal Power Source	Rechargeable Lithium Ion battery, 3.7 V, 2600 mAh	N/A
Charging time for internal lithium ion battery in Charger	Max 3,5 hours, depending on initial state of the battery	N/A
Battery life (fully charged, not connected to mains power)	Min. 3 full charges of 2 hearing instruments, Without hearing instruments: 12 months	N/A
Charging time for Hearing Instrument	Maximum 3 hours, depending on initial state of the battery	Maximum 3 hours, depending on initial state of the battery
Wireless frequency between Hearing Instrument and Charger	2.4 GHz and 333 kHz	2.4 GHz and 333 kHz
ESD tolerance	According IEC 61000-4-2 Electrostatic discharge immunity test standard	According IEC 61000-4-2 Electrostatic discharge immunity test standard
Operating & Charging temperature	0 to 40 °C / 32 to 104 °F	0 to 40 °C / 32 to 104 °F
Storage temperature for charger and Hearing Instrument	-20 to 45 °C / -4 to 113 °F	-20 to 45 °C / -4 to 113 °F

At ReSound we believe life is lived best when it's a natural, effortless flow. This means we take anything that holds you back to heart and will always strive to help you and your clients connect to a life without limitations. We have devoted the past 150 years to advancing lives through technical innovations that, inspired by our Organic Hearing philosophy, place you and your clients at the centre of everything we do, because hearing and hearing care is as individual as a finger-print.

ReSound is part of the GN Group – pioneering great sound from world-leading ReSound hearing aids to Jabra office headsets and sports headphones. The GN Group was founded in 1869, employs over 5,000 people, and is listed on NASDAQ OMX Copenhagen.

Find out how you can help your clients hear like no other with ReSound ONE.

resoundpro.com

-  facebook.com/resoundhearing
-  twitter.com/resoundglobal
-  linkedin.com/company/resound

© 2020 GN Hearing A/S. All rights reserved. ReSound is a trademark of GN Hearing A/S. Apple, the Apple logo, iPhone iPad and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries. Android, Google Play and the Google Play logo are trademarks of Google LLC. Dolby and the double-D symbol are registered trademarks of Dolby Laboratories. The Bluetooth word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.

Worldwide Headquarters

GN Hearing A/S
Lautrupbjerg 7
DK-2750 Ballerup, Denmark
Tel.: +45 4575 1111
resound.com

United Kingdom

GN Hearing UK Ltd.
Kirtlington Business Centre
Portway, Kirtlington
Oxon OX5 3JA
Tel.: +44 1869 352 800
resound.com

Australia

GN Hearing Australia Pty Ltd
Gate C, 19-25 Khartoum Road
Macquarie Technology Park
Macquarie Park NSW 2113
Tel.: (free) 1800 658 955
resound.com

New Zealand

GN Hearing New Zealand Limited
Ground Floor, North Entrance
4 Fred Thomas Drive
Takapuna, Auckland, 0622
Tel.: (free) 0800 900 126
resound.com

Singapore

GN Hearing Pte. Ltd.
456 Alexandra Road,
#22-01,
Singapore 119962
Tel.: +65 6320 9388
resound.com