

ReSound Unite™ Remote Control 2



Product Description

The ReSound Unite™ Remote Control 2 is for clients who appreciate an easy overview of their hearing instrument and accessory settings. They can use it to change programs, adjust the volume, and control streaming from the TV Streamer 2 and Mini Microphone.

The Remote Control 2 allows the client to adjust the volume both monaurally and binaurally. When in streaming mode, they can individually adjust volume of the hearing instrument vs. streamed audio volume.

The Remote Control 2 has a high contrast LCD display so clients can clearly see their instrument settings at a glance.

With the simple keypad, it's easy to choose programs, control volume, and mute hearing instruments. The elegant look and small, handy feel of the Remote Control 2 makes it a pleasure to stay in control.

Easy to connect

- To pair: just press the pairing button on the Remote Control 2 and open and close the battery door of the hearing instruments within 20 seconds
- Pairing can also be done through the ReSound Aventa® 3 fitting software

Easy to use

- At the push of a button you can control volume, mute/unmute, and change programs
- Volume can be adjusted monaurally as well as binaurally
- 'Home' button lets you return to your default settings
- The 'Streaming' button activates streaming from your TV Streamer 2 or Mini Microphone (optional streaming devices)
- Rechargeable with battery status indicator
- Mechanical ON/OFF switch to avoid unintentional activation

2.4 GHz wireless Technology

- A smart, robust, and reliable standard that has a strong, clear, and stable signal

Manufacturer according to FDA:

ReSound North America

8001 Bloomington Freeway
Bloomington, MN 55420
1-888-735-4327
resound.com

ReSound Government Services

8001 Bloomington Freeway
Bloomington, MN 55420
1-800-392-9932
resound.com/veterans

Manufacturer according to
Health Canada:

GN Hearing A/S

Lautrupbjerg 7
DK-2750 Ballerup

Distributor in Canada:

ReSound Canada

303 Supertest Road
Toronto, Ontario
M3J 2M4


1-888-737-6863
resound.com

ReSound

rediscover hearing

Description

- | | |
|---|--|
| 1 ON/OFF switch | 6 Increase volume |
| 2 Display | 7 Decrease volume |
| 3 Streaming (optional TV Streamer or Mini Microphone) | 8 Select left hearing instrument only |
| 4 Home | 9 Select right hearing instrument only |
| 5 Hearing instrument microphone mute | 10 Program change |
| | 11 Pairing |

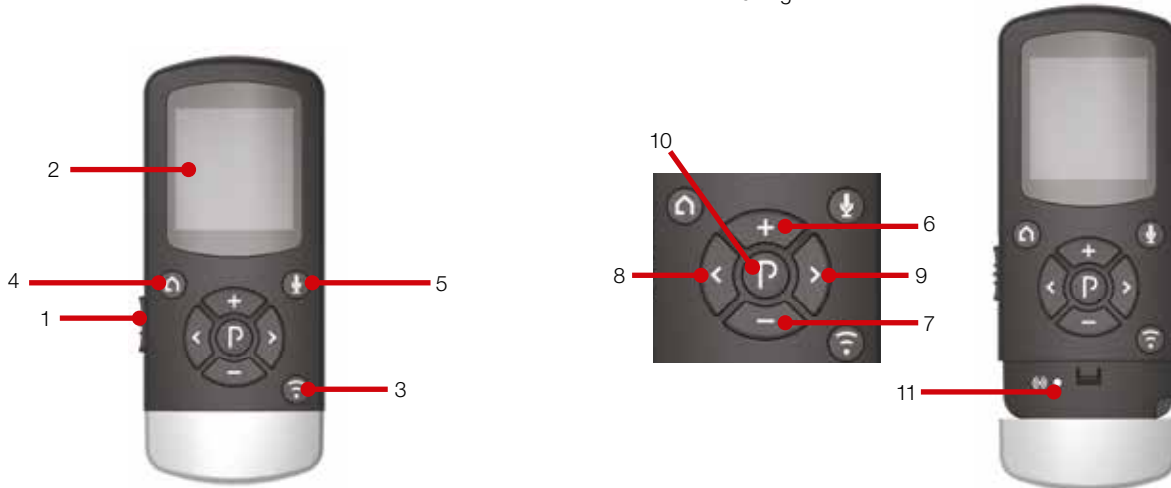

Technical data

| | |
|-----------------------|--|
| Dimensions | 98 x 42 x 15 mm |
| Weight | 35 g |
| Display | 128 x 128 pixel colour display |
| Power supply | Build-in rechargeable battery |
| Battery charge time | Max 2.5 hours using the charger supplied with the unit |
| Wireless connections | 2 hearing instruments per Remote Control 2 |
| Wireless frequency | 2.4 GHz |
| Wireless range | Up to 3 meters |
| Communication speed | 2 Mbit/s |
| Emission | Max. +4 dBm In compliance with all relevant standards |
| ESD tolerance | ±4 kV contact ±8 kV air |
| Operating temperature | 32 to 113 °F |
| Storage temperature | -4 to +113 °F |